

INVACARE®
TDX® SP
POWER WHEELCHAIR



The **Invacare® TDX® SP Power Wheelchair** is designed for people seeking exceptional driving performance, comfort and superior stability in a variety of indoor and outdoor environments.

Call to order

with Invacare®
Power Positioning
1.800.333.6900
invacare.com


with Motion Concepts
Power Positioning
1.888.433.6818
motionconcepts.com



INVACARE CORE TECHNOLOGIES FOR EXCEPTIONAL DRIVING

Center Wheeldrive (CWD)

- Provides the smallest turning diameter and smallest overall footprint vs FWD and RWD for optimum maneuverability
- Provides excellent indoor driving experience due to the exceptional maneuverability

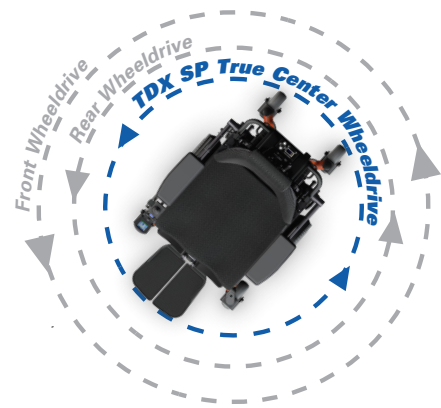

TURNING DIAMETER


SureStep® Suspension and Stability Lock

- Work in synergy to maintain balance and provide stability when transitioning over obstacles
- When maneuvering over everyday obstacles or curbs, SureStep Suspension helps to ensure the seat remains close to level and the power wheelchair maintains virtually constant six wheel contact on the ground
- Stability Lock uses one-way, gas locking cylinders that engage if the rear casters begin to lift off the ground, preventing the chair from lurching forward



G-Trac™ Technology

- Helps provide precise tracking for imprecise environments such as side slopes, uneven or soft ground, or even bunched up carpets or door jams
- Reduces the need for corrective steering actions so the driver uses less effort to reach their destination
- Limiting corrective steering can assist in reduction of fatigue, safer driving, increased driving efficiency, enhanced independence, as well as increased energy to perform other daily living activities

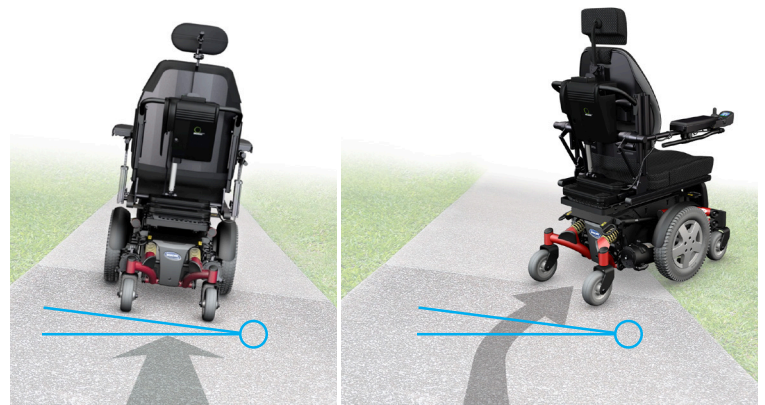

G-Trac: On

G-Trac: Off

MK6i™ Electronics

- Electronics can be programmed to meet individual needs and preferences for powered seating and operation of peripheral devices
- Offers driver-accessible diagnostics capabilities
- Intuitive hand-held programmer offers simple drive programming and setup options



POWER POSITIONING CHOICES TAILORED TO ACCOMMODATE NEEDS FROM SIMPLE TO COMPLEX



ULTRALow
maxx

Premium solution for the most clinically complex needs

Motion Concepts Ultra Low Maxx (ULM) Power Positioning

- Extremely modular system designed to meet an individual's changing needs
- Up to 50° of CG tilt for efficient pressure relief
- Up to 170° of recline with Extended Shear Reduction (ESR) or Power Sliding Back (PSB) to help reduce damaging shear forces and maintain posture
- LNX Power Center Mount Elevating Foot Platform for smooth and effective leg elevation and positioning
- Easy, 'in-clinic' dimensional adjustments to help ensure optimal system set-up and performance
- Seat to floor heights as low as 16.75" with up to 12" of seat elevation to assist with essential, daily activities and promote safe transfers

Motion Concepts' SMART SHIP Program features the ability to order the Invacare TDX SP base, Motion Concepts ULM seating and a full array of ASL options with one order, one invoice.



Full line of Invacare® Matrix® Seating, Cushions and Backs available.

Trusted, value solution for reliability and function


Invacare® Formula™ Power Positioning

- CG tilt allows the seat to slide forward to maintain a constant center of gravity during the tilt cycle
- Rehab (ASBA) seat is adjustable in width and depth
- Up to 166° of recline with Mechanical Shear Management (MSM) maintains consistent positioning and optimally addresses shear forces
- Tilt range 0°-50° or 5°-55°
- Formula systems offer both the Invacare® Contoura™ and Conventional Back styles

Full line of Invacare® Matrix® Seating, Cushions and Backs available with Tilt Only systems when ordered with back post only.

PREMIUM
VALUE

CLASSIC COLORS FOR YOUR PERSONAL STYLE

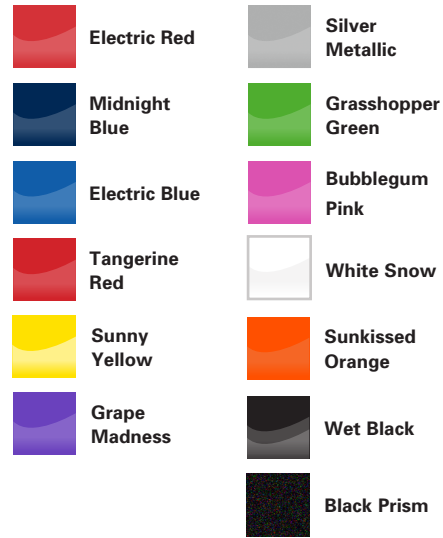

COLORS

HCPCS Code: K0848, K0849, K0856, K0861

Model #	Base Dimensions	Seat Dimensions	Seat-to-Floor Height	Drive Wheels	Ground Clearance	Minimum Base Turning Radius	Product Weight	Weight Capacity	Motor	Max Speed*	Battery	Battery Range*	Limited Warranty**
TDXSP, TDXSP-CG, TDXSP-MCG	35.25" L 24.5" W <i>Narrow</i> 25.5" W <i>Wide</i>	Motion Concepts Ultra Low Maxx 14" W or 18"-22" W 14"-23" D	Invacare Formula or Rehab Seat 16"-22" W 16"-22" D	Motion Concepts Ultra Low Maxx 16.75" 17.25" 18.25" 19.25"	Invacare Formula 18" 19" 20"	Invacare Rehab Seat 16.5" 17.5" 18.5"	Up to 399 lb. <i>Narrow</i> 434 lb. <i>Wide</i>	300 lb.	4-Pole Sealed Housing	Less than 6 mph	Group 24, 22NF	Approximately 18.44 Miles	5 Years Base Frame 3 Years Formula Seat Frame 1 Year Electronics/Motor/ Gearbox

Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Motion Concepts
www.motionconcepts.com

USA
700 Ensminger Rd. Suite 112
Tonawanda, NY 14150
(888) 433-6818
Fax: (888) 433-6834

Canada
84 Citation Dr., Unit 1-7
Concord, ON L4K 3C1
(866) 748-7943
(905) 695-0134
Fax: (905) 695-0138

©2017
Invacare Corporation.
All rights reserved. Trademarks
are identified by the symbols
™ and ®. All trademarks
are owned by or licensed to
Invacare Corporation unless
otherwise noted. Specifications
are subject to change without
notification.
Form No.11-099 Rev. 06/17
170588

