

INVACARE®
ProSPIN™ X4
WHEELCHAIR



The **Invacare ProSPIN X4 Wheelchair** offers the look, feel and performance of rigid wheelchair with the portability and convenience of a folding wheelchair.

Call to order
1.800.333.6900
invacare.com



Features

The uniquely designed single-piece side frame brings aesthetic and performance together in a best possible combination.

The single-piece side frame minimizes both welded and bolt-on connections, strengthening the frame and lessening the amount of propulsion energy lost through flex and vibration. This design also enhances the squareness of the frame.


The curved axle plate completes the flowing side view, making the ProSPIN X4 Wheelchair more like a work of art than a wheelchair.

The durability and ride are legendary, offering all the options a product in this category deserves.

Performance, durability and price puts the Invacare ProSPIN X4 in its own class.

- Fixed front frame version w/ locking flip-up footplate
- Heavy-duty package up to 350 lb.
- Designed from scratch elevating leg rests
- Amputee axle plate option
- Bold black and silver anodized accent package

HCPCS Code: K0005

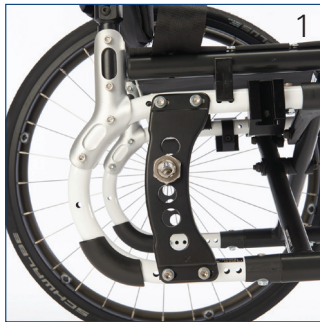

1. Center-of-gravity adjustment
2. Curved handle wheel lock
3. Head tube
4. Swing-in/swing-out front riggings
5. Elevating legrest option
6. Transport Ready Option (TRRO)



Model #	Overall Dimensions	Seat Dimensions	Seat-to-Floor Height	Back Angle	Front Frame Angle	Camber	Tubing Diameter
PROX4S PROX4F70 PROX4F80	Seat width + 8.5" @ 0° camber	14"-20", 22"W 14"-22"D	15.5"- 21.5" Front 13.5"- 20.5" Rear .5" inc	75°-105° 5° inc	60°, 70°	0°, 3°, 6°	1"
	Knee-to-Heel Range	Center-of-Gravity	Wheels	Casters	Total Weight**	Weight Capacity	Limited Warranty
	14"-19.5" 1" inc	Adjustable	20", 22", 24", 25", 26"	3", 4", 5", 6"	30 lb.	250 lb. 350 lb. optional	5 Years Frame/ Crossbraces

Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2018
Invacare Corporation.
All rights reserved.
Trademarks are
identified by the symbols
™ and ®. All trademarks
are owned by or licensed
to Invacare Corporation
unless otherwise
noted. Specifications
are subject to change
without notification.
Form No. 05-005
Rev. 01/18
180124

** 16" x 16", BACK1216, SLIPBACK, 4ALUM2, BRG1, SPCR1, no rear wheels, footrests, anti-tippers, wheel locks or arms.
Please see the Owner's Manual or Invacare Configurator for complete option/specification combinations.

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.